

Bend Coupling

Bend Coupling for Air/Water

A simple solution to quickly connect (and separate) flexible hoses

Available in 2 versions

Version 1 Ø 4,8 mm



Original sizes



Only
1.2 g!

Version 2 Ø 9,0 mm



Version 2 mit Verteiler
with 2, 3, 4 or 5 exits



! All Advantages at a Glance:

- Suitable for connecting flexible hoses
- Service-friendly mounting / dismantling of individual modules and hoses
- Minimal size
- Easy handling (just fix and bend) – no tool required!
- Retrofittable into existing hose systems
- Customizable connection nipples



For further details please see reverse!

All from a single source!

Mytronic has its own development department, manufacture our own tools and moulds and produce over 90 % of our products ourselves. Therefore, we are in a position to respond to any new product requirements quickly and can ensure highest quality standards. Mytronic is certified according to DIN EN ISO 13485.

Mytronic GmbH
 Riedelbach | Technologiepark
 In den Ensterwiesen 7
 61276 Weilrod | Germany
 Tel +49 (0) 6083 95989-0
 Fax+49 (0) 6083 95989-11
 info@mytronic-dental.com

Bend Coupling

Bend Coupling for Air/Water

A simple solution to quickly connect (and separate) flexible hoses


The Mytronic bend coupling is a simple solution to quickly connect and separate flexible hoses.

As the bend coupling (version 1) is made of robust plastic material, it can operate with pressures of up to 8 bar at room temperature.

Functioning of the Bend Coupling

For opening the bend coupling, bend it over the central axis and outer diameter.

For closing, fix it diagonally in slot and pull it up with slight pressure to click into place.

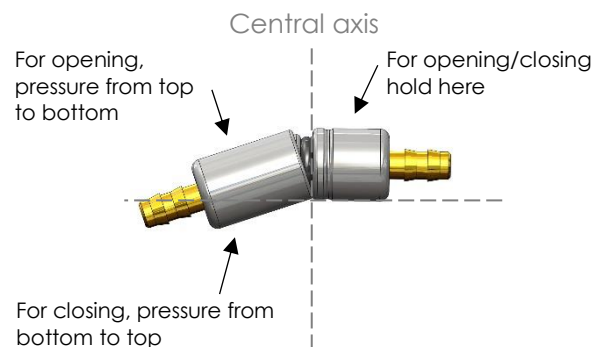

Accessoires

- Safety sleeve for bend coupling version 1



Application example – Connecting individual modules



All from a single source!

Mytronic has its own development department, manufacture our own tools and moulds and produce over 90 % of our products ourselves. Therefore, we are in a position to respond to any new product requirements quickly and can ensure highest quality standards. Mytronic is certified according to DIN EN ISO 13485.

Mytronic GmbH
Riedelbach | Technologiepark
In den Ensterwiesen 7
61276 Weilrod | Germany
Tel +49 (0) 6083 95989-0
Fax+49 (0) 6083 95989-11
info@mytronic-dental.com