


REDUCING AEROSOLS

REGULATION OF WATER AND SPRAY MIST ON MYTRONIC TURBINE HOSES

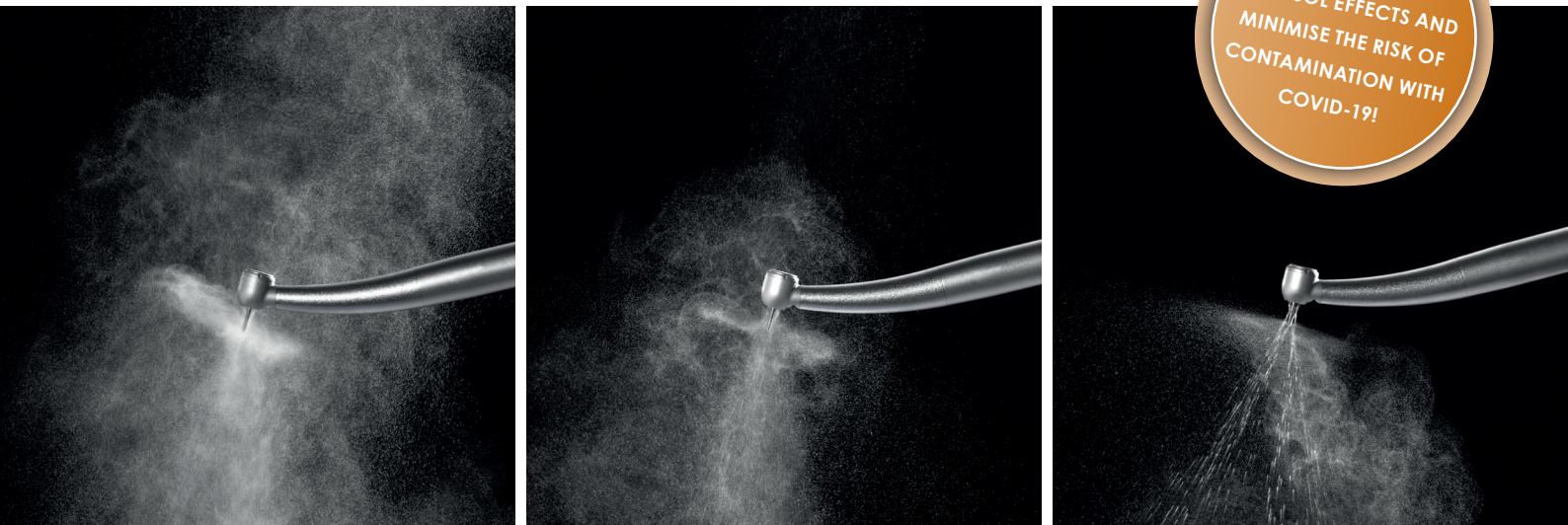

HOW CAN YOU PROTECT YOUR CUSTOMERS AND THEIR PATIENTS AGAINST POTENTIAL DANGER DUE TO AEROSOLS?

Use our tried and tested turbine hoses with the new water and spray mist regulation system.

The goal is to reduce spray mist with this integrated function and to ensure an appropriately adjusted spray supply at the same time, or to just generate a flow of pure water*. The cooling function remains untouched.

Minimise the risk of contamination with the Covid-19 virus by using this special solution from Mytronic.

*(in accordance with ISO 14457:2017-10)

SPECIAL FEATURES

- Reduction of aerosols through water & spray mist regulation directly on the sleeve on the instrument side
- Proportional regulation of water & spray mist with unchanged performance or cooling
- Available for hose connections in accordance with the standard DIN EN ISO 9168:2009-10 (Hose connectors for air driven dental handpieces)
- Simple operation during treatments
- Suitable for use on all commercially available dental units
- Testing option via adapter hose (demo)



For more information please visit:
www.mytronic-dental.com

Would you like to find out more about the water & spray mist regulation system from Mytronic?
 Please call us under **+49 (0) 6083 95989-20**.

Mytronic GmbH

Riedelbach | Technologiepark
 In den Ensterwiesen 7
 61276 Weilrod | Germany

Tel +49 (0) 6083 95989-0
 Fax +49 (0) 6083 95989-11
info@mytronic-dental.com